

Blood Lactate Meter

Lactate Pro 2

LT-1730



For quantitatively measuring the lactate level
15 sec from blood detection

330 test results can be stored up (Up to 3 users)

Blood Lactate Meter

Lactate Pro 2

LT-1730



Small sample volume: 0.3 µL

Only 0.3 µL blood sample is needed for measurement.

Measurement speed:

15 sec from blood detection

Compact - light & handy

Meter dimensions are 100 (L) × 50 (W) × 12 (H) mm.

Meter weight is just 45 g (including battery).

User identification

Up to 3 users can individually save their results up to 330 tests in total.



Sample volume

- 0.3 µL
- 1.0 µL
- 5.0 µL




Test strip for Lactate Pro 2

Test strips are individually foil packed.

Measurement flow



Insert new test strip.

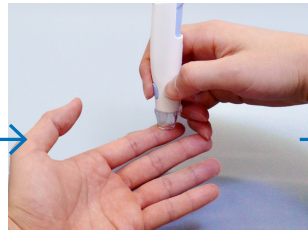

Collect blood.



Apply blood drop to start measurement.



Blood lactate levels are then displayed.

Specification

Test item	Blood lactate level	Battery	3 V lithium battery (CR2032) × 1
Sample	Fresh capillary whole blood	Memory capacity	330 test results (A maximum of 330 test results for all three measurer codes combined.)
Sample size	0.3 µL	Dimensions	50 mm (W) × 100 mm (H) × 12 mm (D)
Test strips	Lactate Pro 2 Test Strip		
Test range	0.5 to 25.0 mmol/L (5 to 225 mg/dL)		
Test time	15 seconds from blood detection		

Note : Not all products shown on this catalogue are offered in every market. Product features may vary by country and are subject to change. For products specific to your area, please visit the website appropriate to your country or contact ARKRAY for more information.

Legal manufacturer

arkray factory, inc.

1480 Koji Konan-cho, Koka, Shiga-ken 520-3306, JAPAN

European representative

arkray europe, b.v.

Prof. J.H. Bavincklaan 2, 1183 AT Amstelveen, THE NETHERLANDS

<https://www.arkray.eu/english/>

*Designs and specifications may be changed without prior notice.

arkray global business, inc.

Yousuien-nai, 59 Gansuin-cho, Kamigyō-ku,

Kyoto 602-0008, JAPAN

TEL +81-75-662-8979 FAX +81-75-431-1202